

### **Cascade cooling system**



Cascade system ULT freezer that consists of two stages cascade cooling technology with absolute hydrocarbon refrigerant, which are R290 and R170 as natural gas. The inverter compressors programmed to help energy savings under long period preservation or provide fast temperature recovery while frequently door open to maintain specimens liveness.



**5.25** 



Energy consumotion

Pulling down to -80°C

Temp. recovery after door opened for 1 min

Features

 $^*$ Typical freezer data based on internal testing with freezer set point at -80°C and ambient temperature at 25°C.

#### **DW-HL778S**



Cascade system
Two stages cascade system



Touch screen 10"LCD touch screen



Hydrocarbon refrigerant Environment-friendly natural gas



Inverter technology
Two inverter compressors



Vaccum insulation
VIP insulation foaming




Alarms system 10 kinds of alarm functions


### Structure



# Dimensions

# Specification

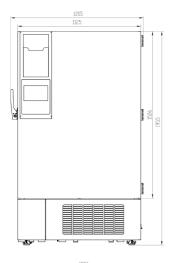

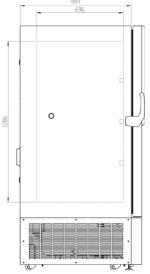

-86°C Ultra-low Temperature Freezer	
Model	DW-HL778S
Capacity	778L/27.5Cu.Ft.
Internal size (W*D*H) mm/inch	865*696*1286/34.0*27.4*50.6
External size (W*D*H) mm/inch	1205*1025*1955/47.4*40.4*77.0
Package size (W*D*H) mm/inch	1320*1155*2171/52.0*45.5*85.5
NW(kgs)	330/420
Performance	
Temperature range	-40∼-86°C
Ambient temperature	16-32°C
Cooling performance	-86°C
Climate class	N
Controller	Microprocessor
Display	10" HD Intelligent Touch Screen
Refrigeration	
Compressors	2pcs
Cooling method	Direct Cooling
Defrost mode	Manual
Refrigerants	High stage: R290, Low stage: R170
Insulation thickness(mm)	130
Construction	
External material	Cold rolled steel sheet
Inner material	Stainless steel
Shelves	3(stainless steel)
Access test port	3pcs. Ø 25 mm
Casters	4+(2 leveling feet)
Data Logging/ Interval /Recording time	USB/Record every 5 minutes / 10 years
	Remote alarm contact, RS485
Interface	Yes
Backup battery Alarms	
Temperature	High the section of t
remperature	High/Low temperature, High ambient temperature
Electrical	Power failure, Low battery, Battery power detection
Liectricat	failure
System	Sensor failure, Main board communication error,
Accessories	Condenser overheating alarm, Door ajar
	714 004 5500 07
Freezer Rack Model	ZKL304-552B-CT
Freezer Rack/Box (Quantity)	24/600
et at teat	
Electrical Data	110 2200//500 50Hz
Electrical Data Power supply(V/HZ) Rated current (A)	110-230V/50&60Hz 14.76A

## Options

#### Cascade models support following optional features:

- Modbus protocol for data transfer
- Card access or finger print access
- O Comform to 21 CFR 11 regulation
- O CO2 back up system

- Inkless temperature chart recorder
- O Racks & Boxes & Tubes

+86 55164317378

www.melingbiomedical.com

overseas.marketing@zkmeiling.com

No.1862 Zishi Road, Hefei, Anhui, 230601, China

